

RELIABLE HEARING IN ANY SITUATION

RELY ON REXTON FOR AN ACCURATE
ALL-ROUND HEARING EXPERIENCE



REXTON

OUR LATEST INNOVATIONS FOR RELIABLE HEARING



BiCore technology




BiCore hearing aids



MotionCore technology



M-Core hearing aids

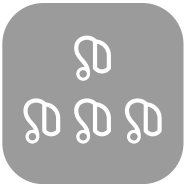

Hearing aid model overview

BICORE C R-Li and C R-Li T 6

The most popular RIC, now with BiCore technology

M-CORE B 15

Rechargeable BTE power, with MotionCore

M-CORE R 20

Our trusted RIC with MotionCore technology

M-CORE SR 22

Award-winning design – with exchangeable receivers

M-CORE iX 23

Our smallest instant-fit CIC, with MotionCore

Adore Custom Products 24

A complete line of custom products

Model overview 26

Proven technology



BICORE TECHNOLOGY



Having a conversation in busy and noisy environments is a challenge for hearing aid users.

They need to focus attention on speech yet still stay aware of the sounds around them. Not being able to do so can be confusing, tiring and even dangerous. Our Rexton BiCore hearing aids solve this challenge. BiCore technology automatically splits the focus area from the background sound and then processes the two signals separately. When the signals are merged back together the result is a smooth and fuller listening experience for the users even in demanding environments.



REXTON APP: A HEARING AID WEARER'S ALL-IN-ONE SMARTPHONE COMPANION

The Rexton App allows wearers to control their hearing aids conveniently and discreetly. It may also be used for remote tuning adjustments. The app also features our new Rexton Assist, for self-service hearing aids programming.

Download the Rexton App now!



Rexton Assist

Rexton Assist, accessed via the Rexton App, works on a live, learning network that continuously evolves to provide solutions to any challenge faced by the hearing aid wearer.

Its initial training was performed using 13 million data points, learning from every interaction and getting better and more accurate all the time.

Rexton Assist brings convenience for both the user and the dispenser, improving the hearing aid performance based on the needs of the wearer, at any time, in any place.



Made for iPhone and Android (ASHA)

BiCore hearing aids can stream calls and audio from iPhones and Android phones.*



Streaming

You may direct stream from your TV and other audio devices using the Smart Transmitter 2.4 or the Smart Mic.

*Android phones must operate on ASHA protocol.



BICORE C R-LI AND C R-LI T

OUR MOST RELIABLE HEARING SOLUTION YET
WITH THE LATEST BICORE TECHNOLOGY

The new BiCore C R-Li and C R-Li T are rechargeable and also feature Rexton Assist, combined with the 3rd generation of motion sensors and with My Voice technology.

Available colors



BICORE C R-LI and C R-Li T







Performance Level	80
Battery	Lithium-ion battery
	Charging time 4 h
	BICORE C R-LI Battery runtime (without streaming) up to 28 h Battery runtime (incl. 5 h streaming) up to 24 h
	BICORE C R-LI T Battery runtime (without streaming) up to 39 h Battery runtime (incl. 5 h streaming) up to 36 h
FOG in dB	46/60/65/75
Accessories	Standard or Charging+ Station
	Smart Mic
	Smart Transmitter 2.4
	Smart Key
	Rexton App

Features	BiCore Processing
	Rexton Assist
	Motion Sensor
	IP 68 Rated
	Soundpro
	My Voice (own voice processing)
	Direct Streaming
	Auto Volume
	Wireless Sync
	Directionality
	Noise Reduction
	Wind Noise Reduction
	Reverb Reducer
	Music Enhancer (Live/Recorded/Playing)
	T-Coil (BiCore C R-Li T)
	XPhone
Programming	NoahLink Wireless

Chargeable via



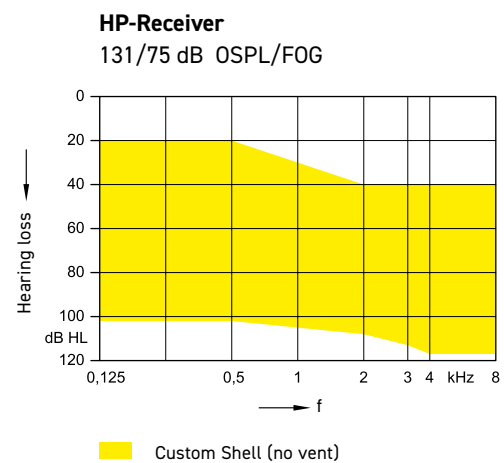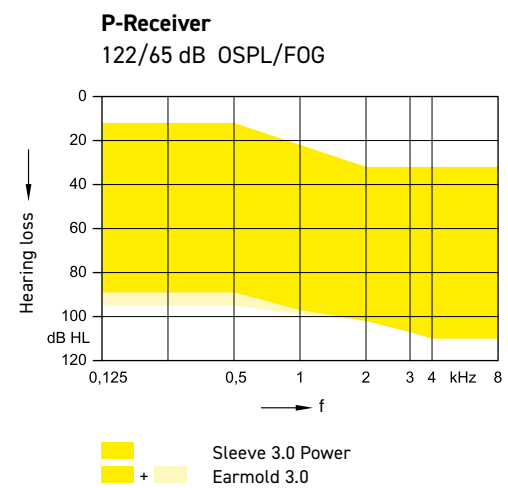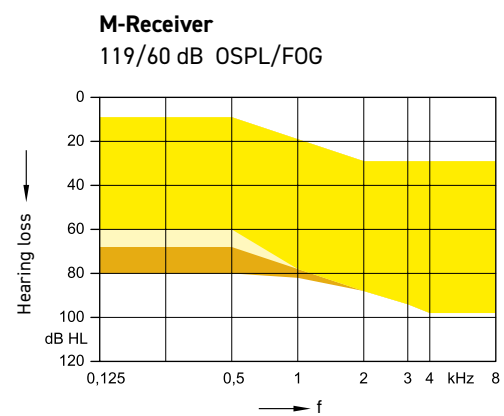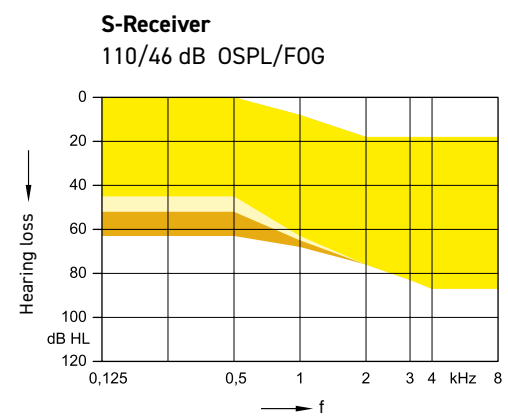

CHARGERS

Connexx Travel Charger

This portable charger provides up to 3 charges on the go – without connecting to a power outlet.

- Charges in 4 hours
- 30-minute fast charging
- Easy insertion/removal of aids
- Automatic on/off
- USB-C, power adapter

Available for

- Standard: BiCore C R-Li and C R-Li T
- Standard: M-Core SR



Connexx Charging+ Station

A multi-function charger that makes taking care of your hearing aids easy. It not only charges your hearing aids but also includes a drying function to ensure there is no moisture build-up. In addition there is a UV sanitizing cycle that helps destroy any harmful germs or bacteria from the surface areas.

- Charges in 4 hours
- 30-minute fast charge
- Integrated drying cycle
- Integrated UV sanitizing cycle

Optional for

- BiCore C R-Li/Li T
- M-Core B Li P/HP



Connexx Standard Charger

The stationary Standard charger is a very stable and easy to use unit, with a robust design and construction.

- Charges in 4 hours
- 30-minute fast charging
- Easy insertion/removal of aids
- Automatic on/off
- Fixed charging cable
- Built in magnet holds hearing aids in place

Available for

- BiCore C R-Li/Li T (optional)
- M-Core B-Li P/HP (included)



Connexx Charging Station

- Charges in 4 hours
- 30-minute fast charging
- Integrated drying cycle

Available for

- M Core B Li M (included)



EARWEAR 3.0 FOR BICORE RICS



MINIRECEIVER 3.0

Easy to insert and remove. No tool necessary.

Lengths	0, 1, 2, 3, 4
Power level	S, M, P, HP
Concha Lock	One size fits all




MEASURING TOOL 3.0

Quickly and accurately determine the best length for miniReceiver.



EAR TIPS 3.0

Soft tips for a comfortable fit.

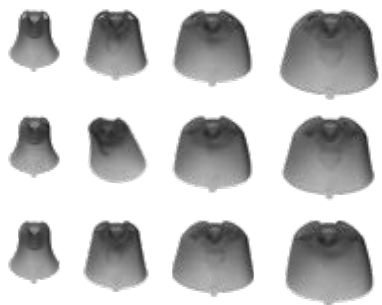

DOMES	
Open	4mm, 7mm, 10mm
Tulip	8mm, 12mm
SLEEVES*	
Open	XS, S, M, L
Closed	XS, S, M, L
Power	XS, S, M, L

* Coming May 2021



NANOCARE 3.0 WAXGUARDS

Keep wax from damaging receivers with replaceable waxguards.

MOTIONCORE TECHNOLOGY



Most hearing aids analyze hearing situations solely on acoustic criteria, then classify them into a small range of predefined acoustic environments. But when people rely on their senses to keep them safe and sound, hearing clearly in every situation is a must – especially when they're on the job and moving through the day, where environments can change quickly. Our M-Core hearing aids use proven MotionCore technology that automatically recognizes and adapts to different hearing situations. This ensures users can hear speech clearly and get the most reliable impression of their environment.

UP TO 20 TIMES MORE HEARING SITUATIONS THAN BEFORE, AUTOMATICALLY



High Resolution Soundpro

Using a combination of inputs including real-time motion data, signal to noise ratio, and sound directionality, High Resolution Soundpro can continuously adjust the sound settings to provide the most reliable impression of the user's environment, automatically.



Motion Sensor

MotionCore technology features an integrated real-time motion sensor that detects a user's movement, and uses the data to even more accurately determine their sound situation. With this, M-Core hearing aids can automatically detect up to 127 different sound situations.



My Voice

The hearing aid recognizes the user's specific voice pattern, and integrates it naturally into what they hear. In turn, this provides the user with a natural-sounding voice that leads to high acceptance rates. And for that the user will communicate more and be more socially active, greatly improving their quality of life.



Telecare


Full remote fitting sessions via Connexx resulting in higher member convenience and satisfaction.

RECHARGEABLE BTE POWER

WITH MOTIONCORE TECHNOLOGY

Our BTE hearing aids come completely rechargeable, even up to our High Power models which have an unbelievable runtime of up to 61 hours on a single charge. They also feature MotionCore technology, which automatically recognizes and adapts to more hearing situations than ever. With our rechargeable high power BTE, more users can automatically get the most reliable impression of their environment than ever. M-Core B-Li hearing aids are designed to be durable, practical and easy-to-use, so your customers can rely on Rexton in any situation.

Available colors for all M-Core BTEs and M-Core RIC 312

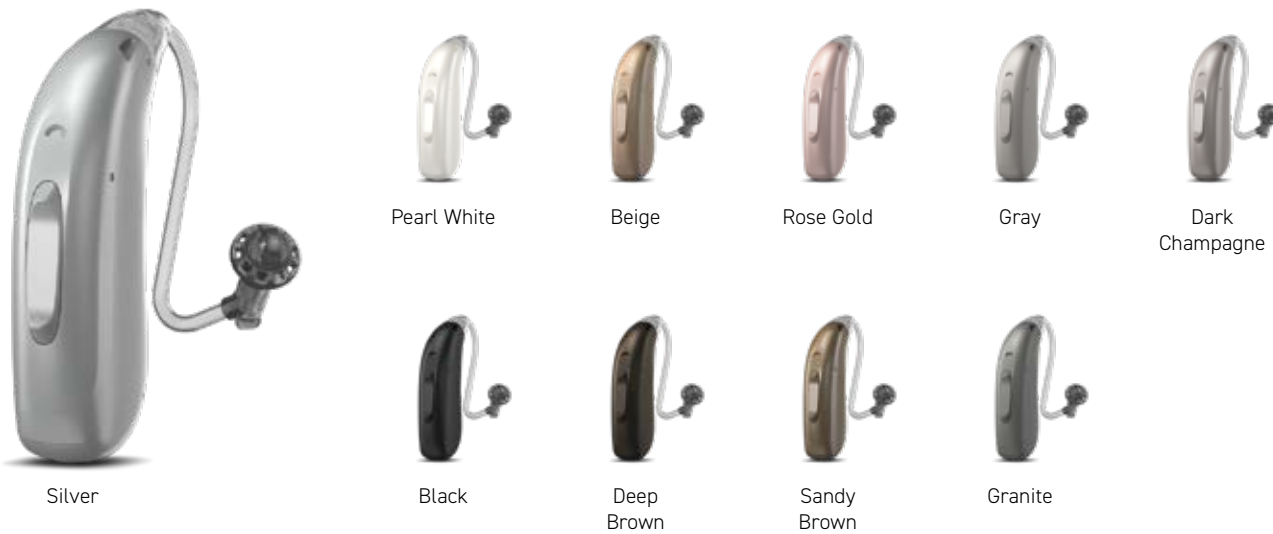


M-CORE B-Li M

Performance Level	80
Battery	Lithium-ion battery
FOG in dB	60
Accessories	Connexx Charging Station
	Smart Mic
	Smart Transmitter 2.4
	Smart Key
	Rexton App
	EarWear 3.0 Thin Tubes and Tips
	Color Conversion Kits
Features	Integrated motion sensor
	Hi-Res Soundpro
	Up to 120 hearing situations
	Direct Audio Streaming
	IP 68 Rated
Programming	NoahLink Wireless

With a compact size that is comparable to a 312 battery device, the M-CORE B-Li M combines a modern design with rechargeability and Bluetooth connectivity.

133/60 dB OSPL/FOG (Earhook)

125/60 dB OSPL/FOG (Thin Tube)

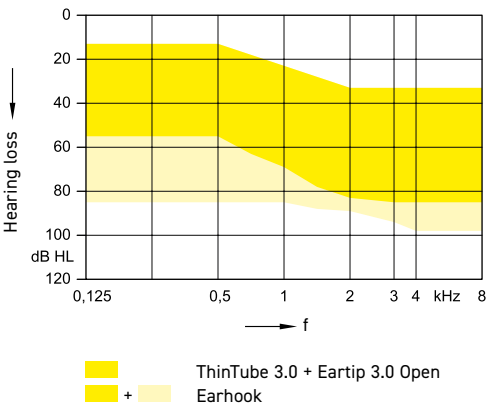

Chargeable via



Connexx Charging Station


RECHARGEABLE BTES

WITH THE POWER OF MOTIONCORE

With a smaller size, and rechargeable battery, the M-CORE B-Li P remains a powerful device that also offers the latest MotionCore technology.

135/77 dB OSPL/FOG (Earhook)

130/66 dB OSPL/FOG (Thin Tube)

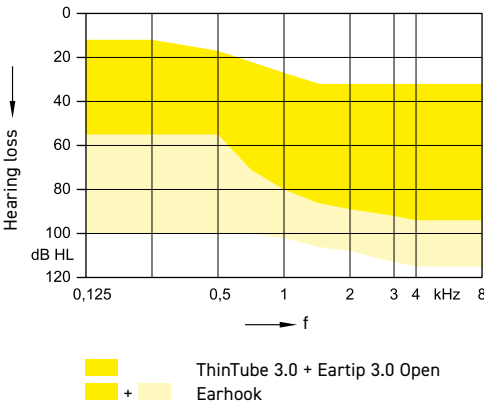

M-CORE B-Li P



Performance Level	80
Battery	Lithium-ion battery
FOG in dB	77
Accessories	Standard or Charging+ Station
	Smart Mic
	Smart Transmitter 2.4
	Smart Key
	Rexton App
	EarWear 3.0 Thin Tubes and Tips
	Color Conversion Kits
Features	Integrated motion sensor
	HiRes SoundPro
	Up to 120 hearing situations
	Direct Audio Streaming
	TeleCoil
	IP 68 Rated
Programming	NoahLink Wireless

Chargeable via



Connexx Charging+ Station
Optional




Connexx Standard Charger

OUR FIRST
RECHARGEABLE
HIGH POWER
BTE

UP TO
61 HOURS
RUNTIME ON
A SINGLE
CHARGE

Our most powerful device, the M-CORE B-Li HP is the smart choice for those with severe to profound hearing loss, combining the latest MotionCore technology with rechargeability.

140/82 dB OSPL/FOG (Earhook)

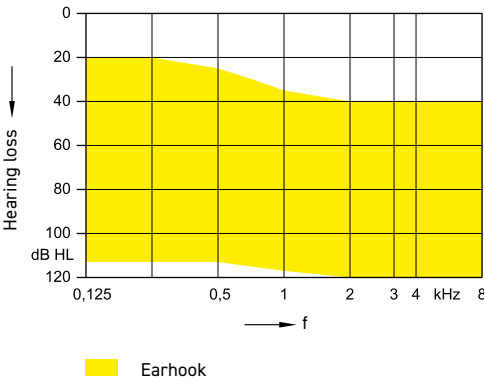

M-CORE B-Li HP



Performance Level	80
Battery	Lithium-ion battery
FOG in dB	82
Accessories	Standard or Charging+ Station
	Smart Mic
	Smart Transmitter 2.4
	Rexton App
	Color Conversion Kits
Features	Integrated motion sensor
	HiRes SoundPro
	Up to 120 hearing situations
	Direct Audio Streaming
	TeleCoil
	IP 68 Rated
Programming	NoahLink Wireless

Chargeable via



Connexx Charging+ Station
Optional



Connexx Standard Charger

EARWEAR 3.0 FOR BTES

M-Core BTEs have new accessories to help personalize them for the wearer when wanting to use thin tubes instead of regular tubing and earmolds.

- Available for M-Core B-Li M and M-Core B-Li P

THIN TUBE 3.0



Lengths	0, 1, 2, 3, 4
Diameter	0.9 (Standard) 1.4 (Power)
Concha Lock 3.0	One size fits all

MEASURING TOOL 3.0

Quickly and accurately determine the best length for miniReceiver.

EAR TIPS 3.0

Soft tips for a comfortable fit.

DOMES	
Open	4mm, 7mm, 10mm
Tulip	8mm, 12mm
SLEEVES*	
Open	XS, S, M, L
Closed	XS, S, M, L
Power	XS, S, M, L

* Coming May 2021


NANOCARE 3.0 WAXGUARDS

Keep wax from damaging receivers with replaceable waxguards.



OUR TRUSTED RIC

WITH MOTIONCORE TECHNOLOGY

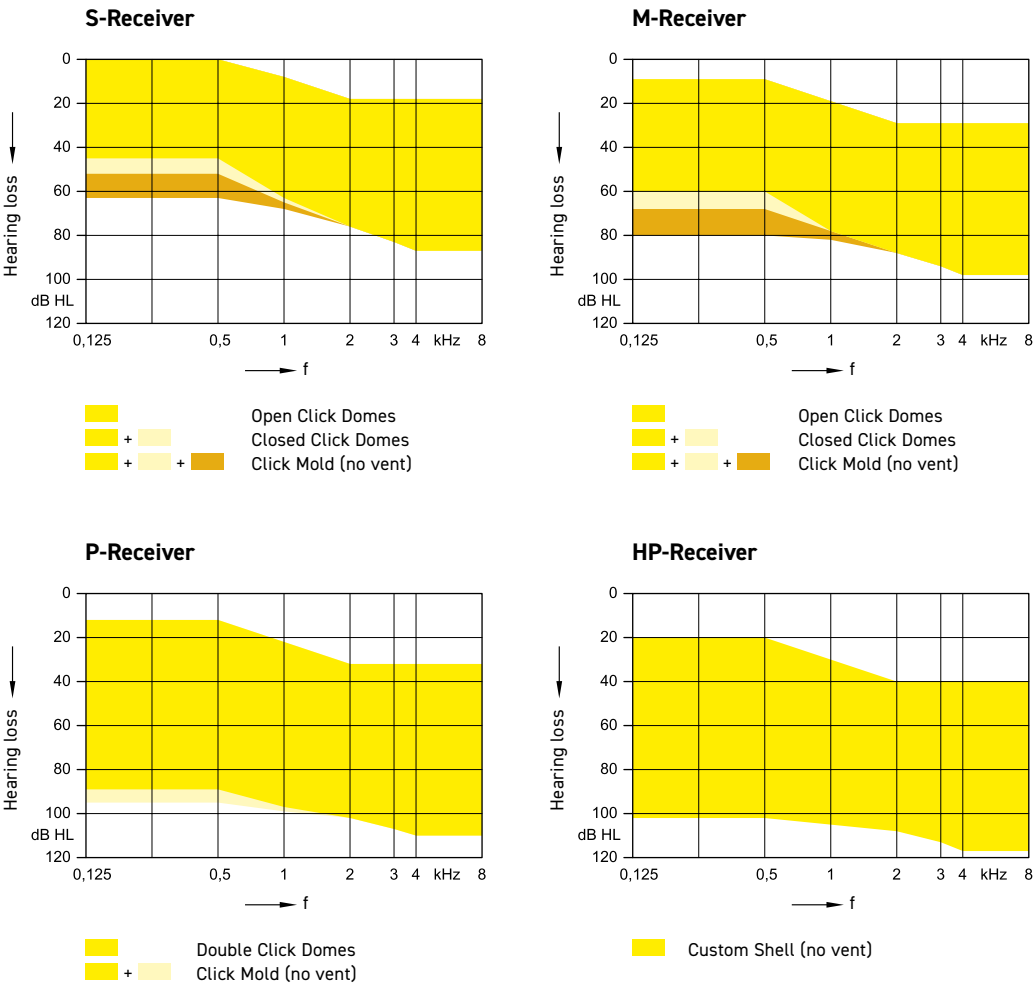

M-CORE R312



Performance Level	80
Battery	312
Accessories	Smart Mic
	Smart Key
	Smart Transmitter 2.4
	Rexton App
Features	Integrated Motion Sensor
	Programmable rocker switch
	Color Conversion Kits
	Color Conversion Kits with design integrated telephone coil (M-Core R312T)
	IP 68 Rated
Programming	NoahLink Wireless

M-Core R312 is available with an integrated Telecoil for optimal hearing where loop systems are available.






M-Core R-312T with optional telecoil

AWARD-WINNING DESIGN

WITH EXCHANGEABLE RECEIVERS

Our award-winning Slim-RIC design features proven MotionCore technology that automatically recognizes and adapts to up to 20 times more hearing situations than before.

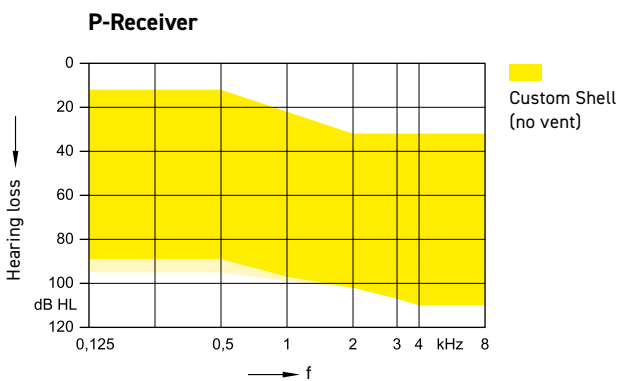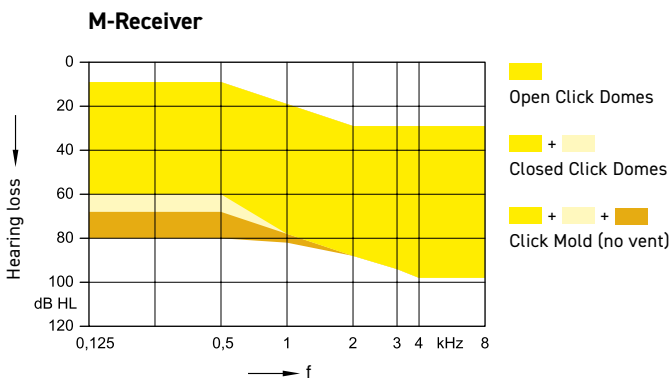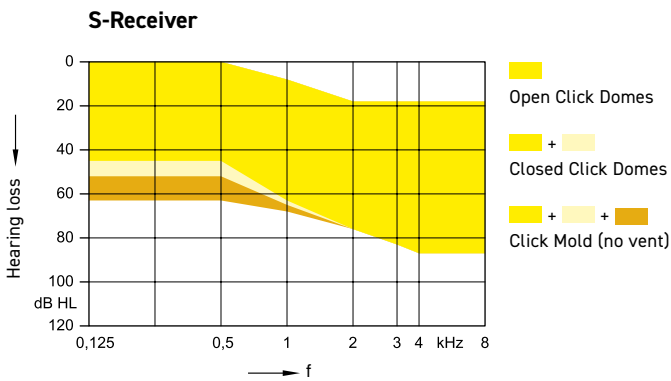

M-CORE SR


Performance Level	80
Battery	Lithium-ion battery
Accessories	Connexx portable charging case
	Smart Mic
	Smart Transmitter 2.4
	Smart Key
	Rexton App
Features	Integrated Motion Sensor
	127 sound situations
	Exchangeable Receivers
	Direct Streaming
	IP 68 Rated
Programming	NoahLink Wireless

Chargeable via



Connexx Travel Charger

Choose from 5 distinct color combinations

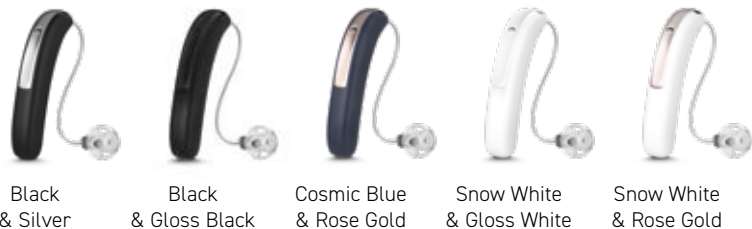


INSTANT FIT

WITH THE HIDDEN POWER OF MOTIONCORE



Our virtually invisible CIC design features proven MotionCore technology that automatically recognizes and adapts to up to 10 times more hearing situations than before. This ensures users hear speech clearer and get the most reliable impression of their environment, automatically. M-Core iX hearing aids are designed to be discreet, practical and easy-to-use.

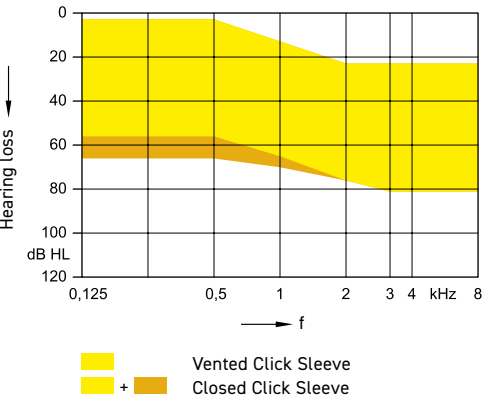


M-CORE iX

Performance Level	80
Battery	10
Accessories	Smart Key
	Rexton App
Features	HD Soundpro
	70 sound situations
	Instant-Fit Sleeves
Programming	Flex connector HiPro, Connexx Air

Easy to fit – no impressions required
Use with all Connexx domes and sleeves




Accessories

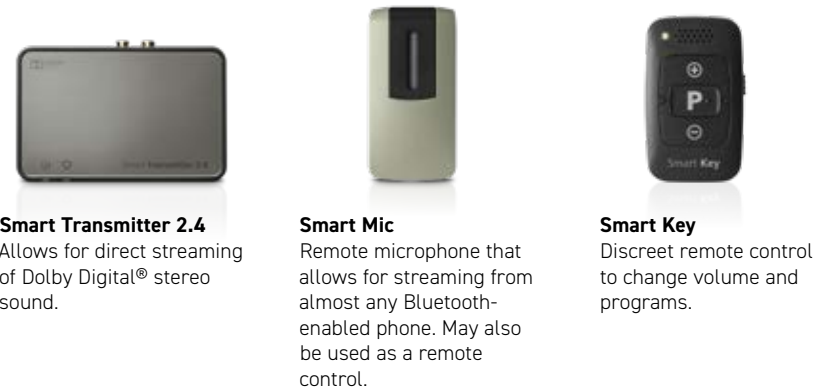

ULTRA-DISCREET WITH MYCORE FEATURES



Our smallest custom models offer customers and discretion and performance with our MyCore processing featuring Single Mic Directionality.

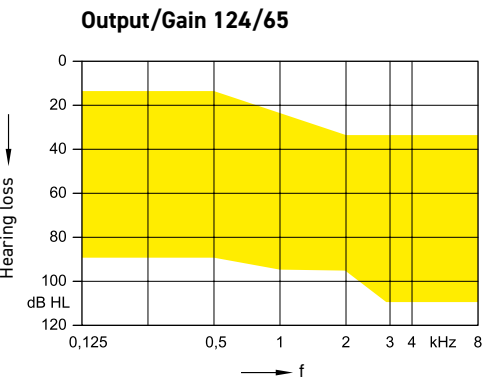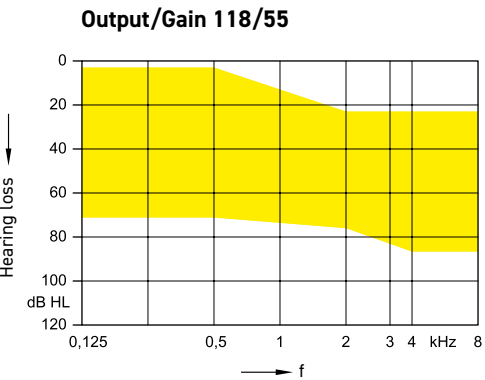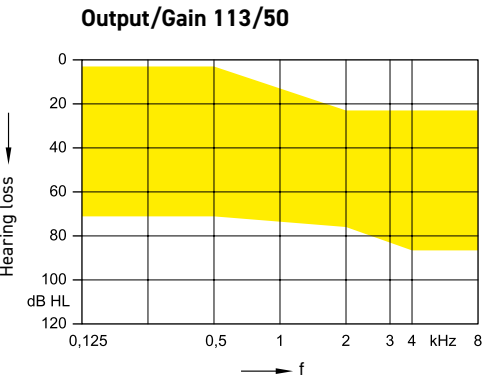


ADORE IIC & CIC






Performance Level	80
Battery size	10
FOG/OSPL in dB	113/50 (CIC and IIC) 118/55 124/64 (CIC only)
Accessories	Smart Key Rexton App
Features	SoundPro Voice Ranger Wireless sync Automatic/Adaptive bin. Single-mic Directionality Noise Management SoundSmoothing Directional Wind Noise Reduction Reverb Reducer Bandwidth: Extension/Compression Music Enhancer: Live/Recording Xphone Acclimatization/Data Logging
Programming	ConnexxAir HiPro (CIC only)

Available colors

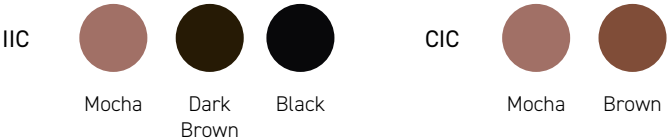

Our full shell and canal models provide excellent speech clarity and noise reduction due to our binaural Wireless Sync. They also offer direct streaming to iOS devices and high power to fit mild to severe losses.

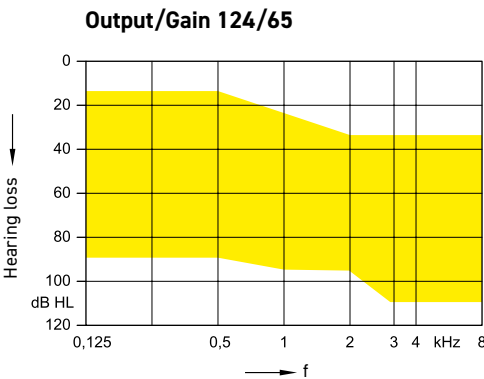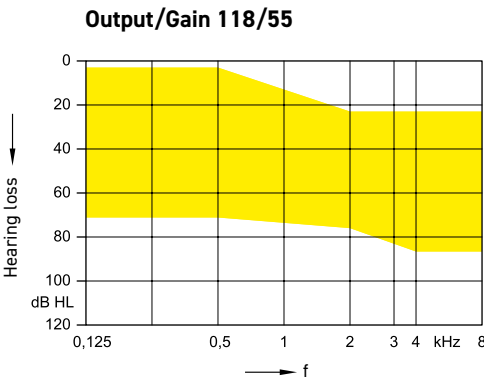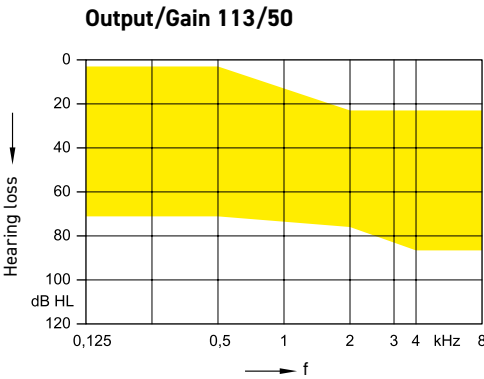

ADORE ITE/ITC








Performance Level	80/80 BT
Battery size	312
FOG/OSPL in dB	50/113 55/118 65/124
Accessories	Smart Mic Smart Transmitter 2.4 Smart Key Rexton App
Features	SoundPro Voice Ranger Wireless sync Directionality Automatic/Adaptive Front & Back, Left & Right, Narrow - Stereo iLock Noise Management SoundSmoothing Directional Reverb Reducer Bandwidth: Extension/Compression Music Enhancer: Live/Recording Xphone Acclimatization/Data Logging
Programming	NoahLink Wireless (BT) ConnexxAir (Non-BT) HiPro

Available colors





A FULL RANGE OF RELIABLE HEARING OPTIONS FOR YOUR CUSTOMERS

	RIC				Slim-Ric			BTE			Instant Fit	Custom			
															
Product	BiCore R-Li NEW	BiCore R-Li T NEW	M-Core R312	Styleline M-Core SR	M-Core B-Li M	M-Core B-Li P	M-Core B-Li HP	M-Core B-Li M	M-Core B-Li P	M-Core B-Li HP	M-Core iX-CIC	Adore IIC	Adore CIC	Adore ITC Bluetooth	Adore ITE Bluetooth
Battery	Li-Ion	Li-Ion	312	Li-Ion	Li-Ion	Li-Ion	Li-Ion	Li-Ion	Li-Ion	Li-Ion	10	10	10	312	312
Programming Device	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	NoahLink Wireless	Connexx Air High-Pro	Connexx Air	Connexx Air High-Pro	NoahLink Wireless High-Pro	NoahLink Wireless High-Pro
TeleCare	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Remote Session	Portal	Portal	Portal	Remote Session	Remote Session
Bluetooth Capable	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	Yes	No	No	No	Yes	Yes
T-Coil Option	No	Yes	Yes	No	No	Yes	Yes	No	Yes	Yes	No	No	No	No	No
Travel Charger	Included	Included	—	Included	No	No	No	No	No	No	No	No	No	No	No
Charging+ Station	Optional	Optional	—	No	No	Optional	Optional	No	Optional	Optional	No	No	No	No	No
Standard Charger	Optional	Optional	—	No	No	Included	Included	No	Included	Included	No	No	No	No	No
Charging Station	No	No	—	No	Included	No	No	Included	No	No	No	No	No	No	No



RELY ON REXTON

It's a tough world out there, and even tougher when your hearing starts to go. At Rexton, we won't let hearing loss hold you back. Good hearing is critical: From getting the job done, to getting home safely and being there for the people who count on you, we know what's at stake. That's why we work hard to deliver proven hearing technology in the most reliable way. We understand what reliability means in real life, and make practical, easy-to-use products you can count on. Our hearing aids have been getting the job done since 1955, so whatever life throws at you, you can rely on Rexton.

www.costco.com/rexton.html

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Sivantos GmbH is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Google Play is a trademark of Google Inc. "Made for iPhone" and "Made for iPad" mean that an electronic accessory has been designed to connect specifically to iPhone or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone or iPad may affect wireless performance. iPad, iPhone are trademarks of Apple Inc., registered in the U.S. and other countries.

Order No. 04537-99C1-7600 © 2020 Sivantos GmbH