

# RELIABLE HEARING IN ANY SITUATION

REACH PRODUCT CATALOG 2026



**REXTON**

# RELY ON REXTON

Rexton hearing aids get the job done every single time. They will perform reliably at home, at work, during exercise, and you can wear them in all weather. That's because they are built tough, using the highest quality materials, and tested rigorously to make sure they can answer all the challenges that life throws at you.

Noisy rooms? No problem, you can focus on what's being said. Drop them accidentally? They are impact proof. Caught in the rain? They are water resistant. Working in a harsh environment? Don't worry about it. And they can handle much, much more. Simply put, they'll keep you prepared for every situation. They are made to be Lifeproof.

RELY ON  
**LIFEPROOF**  
REXTON

# A ROBUST PORTFOLIO FOR RELIABLE HEARING



## Reach Technology

Reach Plus BC .....	6
Reach inoX-CIC Li .....	8
Reach B-Li Rugged .....	10

Chargers .....	12
----------------	----

The comprehensive Rexton App .....	14
------------------------------------	----

Accessories & Apps .....	16
--------------------------	----

EarWear 3.0 & ThinTube 3.0 .....	18
----------------------------------	----

# REACH TECHNOLOGY



Every voice matters, and Rexton Reach makes sure every voice is heard. These rechargeable hearing aids with Multi-Voice Focus enable wearers to keep track of multiple voices simultaneously, even when there is background sound.

## MULTI-VOICE FOCUS

Multi-Voice Focus scans the surrounding environment 1000 times a second with 4 independent focus beams. This enables it to constantly adapt to the changing dynamics of the conversation and make the experience fluid and comfortable. It will boost the prominence of those in the group who are speaking while keeping track of those who are silent, and all without switching between modes, and with no delays.

## VOICE STABILIZER

Voice Stabilizer controls the compression to ensure all voices within the focus beam are amplified in the best linear ways when the sound gets too loud. It allows for more aggressive compression of the noise and more linear processing of speech to ensure a strong contrast between speech and noise. Benefit: voices remain untouched while surrounding noise is managed.

## CONNECTIVITY FOR EVERYONE WITH BLUETOOTH CLASSIC

Reach Plus BC hearing aids connect to 99% of devices via Bluetooth® Classic, for simple, reliable pairing and streaming connections. Reach Plus BC is also Auracast Ready with a future over-the-air Firmware update.

## LIFEPROOF

Rexton hearing aids are built using robust components and are tested extensively to ensure they can endure harsh environments and you can rely on them in everyday challenges.



## SOUNDPRO 3.0

SoundPro 3.0 automatically steers all sound processing features, including the new Multi-Voice Focus\* inputs, based on the analysis of the acoustic scene. When changing the strategy to Natural Focus or Speech Focus, the recommended balance of features is modified. Balanced Focus is recommended for enhanced automatic performance.

\* Only for 2-Mic devices and binaural fittings



## BLUETOOTH AUDIO STREAMING

Reach hearing aids can stream calls and audio from virtually any kind of Bluetooth-enabled smartphone.\*\*



## TELECARE

Full remote fitting sessions via Connexx resulting in higher wearer convenience and satisfaction.



## HEADSET MODE

Rexton Reach can use the hearing aid microphones as a headset so wearers can enjoy hands-free phone calls.\*\*



## MY VOICE 2.0

Automatically reduces the sound of the wearer's voice when they speak so their voice sounds more natural to them.



## AUTOECHO REDUCER

Automated environment scan to reduce effects of disturbing reverberation on the perceived sound signal.



## ATEX

Rexton Reach have ATEX certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments.\*\*\*

\*\* Reach Plus BC: 99% device compatibility as wireless streaming is via Bluetooth Classic. BT streaming and headset mode not available for Reach inoX.

\*\*\*They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills. The ATEX group/class is printed on the ATEX card, which must be available at hand when the hearing aids are used in explosive atmospheres. Ask your hearing care professional for your ATEX card.



# REACH PLUS BC

## REACH OFFERS A SOLUTION FOR EVERYONE WITH BLUETOOTH CLASSIC



Reach R-Li T BC

- › Enabling users to be fully involved in group conversations with Multi-Voice Focus
- › Handsfree phone calls via Headset Mode with Call Control when connected via Bluetooth Classic
- › Portable charger with built-in power bank
- › Lifeproof to withstand everyday challenges
- › ATEX certified
- › AI with Rexton Assist helps you to seamlessly customize your hearing aids

### Chargeable via




Travel Charger RIC Included



Charging+ Station R Optional



Fitting ranges are based on 2 ccm coupler:

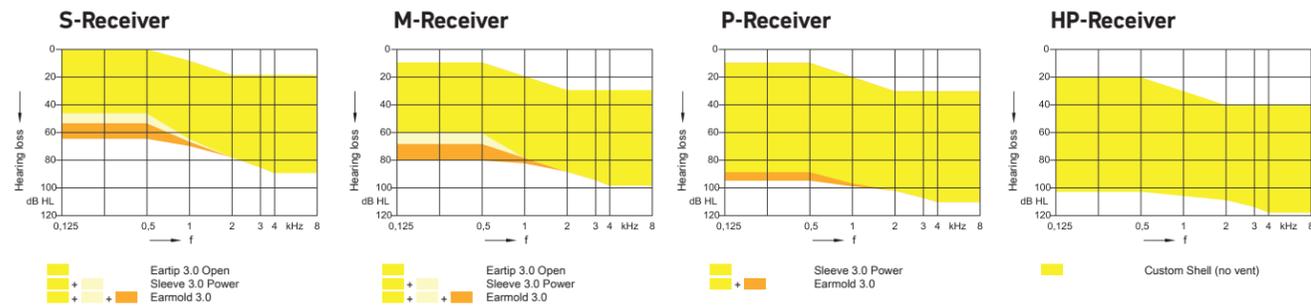

- 46 dB / 110 dB SPL (2 ccm coupler)
- 56 dB / 120 dB SPL (Ear simulator)
- 60 dB / 119 dB SPL (2 ccm coupler)
- 70 dB / 129 dB SPL (Ear simulator)
- 65 dB / 122 dB SPL (2 ccm coupler)
- 75 dB / 131 dB SPL (Ear simulator)
- 75 dB / 131 dB SPL (2 ccm coupler)
- 83 dB / 138 dB SPL (Ear simulator)

## REACH PLUS BC

Features / Performance Level	80
Channels/Controls/Programs	48/20/6
LED	✓
T-coil	✓
Bluetooth	✓
Wireless sync	✓
IP68	✓
XPhone	-

Battery	
Lithium-ion battery	✓
Charging time	4 h
Battery runtime (without streaming)	up to 54 h
Battery runtime (incl. 5 h streaming)	up to 36 h

Technology Highlights	
Multi-Voice Focus	✓
Focus Beams	4
Voice Stabilizer	✓
Speech Preservation Technology	✓
My Voice 2.0	✓
Motion sensor for scene analysis	✓

Directionality	
Automatic adaptive	✓
iOmni Sound Locator	✓
Auto iFocus 360	✓
iFocus 360	✓
Auto Stereo iLock	✓

Noise Reduction	
Noise Management	✓
Sound Smoothing	✓
Wind Noise Reduction	✓
AutoEcho Reducer	✓
Reverb Reducer	✓

Sound Quality & Control	
Bandwidth extension	✓
Bandwidth compression	✓
Intelligent Feedback Preventer	✓
Music Enhancer (presets)	3
Automatic Acclimatization	✓

Audio Streaming	
Bluetooth Classic	✓
LE Audio	-
MFi	-
ASHA	-

Headset Mode	
Bluetooth Classic	✓
LE Audio	-
MFi	-

TeleCare	
Full live remote tuning	✓
Basic portal remote tuning	✓

Accessories & Apps	
Smart Key	✓
TV Sound	✓
Rexton App	✓

Chargers	
Travel Charger RIC	✓
Charging Station R	✓
Charging+ Station R	✓

Programming	
Noahlink Wireless	✓
ConnexAir	-
HI-PRO	-

✓ available - not available



### ATEX CERTIFICATION

Reach hearing aids have an ATEX certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments. They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills.

# REACH INOX-CIC LI

The first invisible hearing aid with rechargeability that is compact, and instantly fits the user's ear canal comfortably. Reach inoX-CIC Li allows the wearer to be fully involved in the conversation without others paying attention to their hearing aid.



- > Virtually invisible design for utmost discretion
- > Rechargeable Li-ion battery for full day use (~28H)
- > Rexton App compatible
- > 4 charges on-the-go with built-in power bank
- > Programmable via ConnexxAir
- > ATEX certified
- > Water resistant, IP68
- > Lifeproof
- > Compatible with Instant Fit CIC Earmold 3.0

## Faceplate Colors




## Chargeable via



IF CIC Travel Charger Included

All fitting ranges are based on 2 ccm coupler:

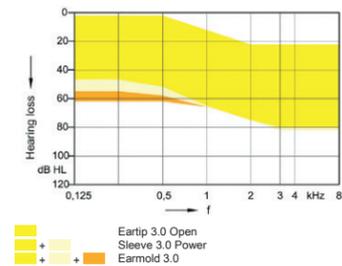

- 50 dB / 114 dB SPL (2 ccm coupler)
- 61 dB / 126 dB SPL (Ear simulator)

## REACH INOX- CIC LI

<b>Performance Level</b>	80
Channels/Controls/Programs	48/20/6
<b>Battery</b>	
Lithium-ion battery	✓
Charging time	4 h
Wearing time	up to 28 h
<b>Matrix</b>	114/50
<b>Channels</b>	48
<b>Bands</b>	20 Adjustable
<b>Programs</b>	6
<b>Accessories and Apps</b>	
Smart Key	✓
Rexton App	✓
<b>EarWear 3.0</b>	
Domes & Sleeves	✓
Instant Fit CIC Earmold 3.0	✓
<b>Charger</b>	
IF CIC Travel Charger	✓
Qi Charging	✓
Holds up to 4 charges	✓
<b>Programmer</b>	
ConnexxAir	✓

Features	
Soundpro 3.0	✓
Voice Stabilizer	✓
OneMic Directionality	✓
Wireless Sync 4.0 (4x speed of BiCore)	✓
AutoEcho Reducer	✓
Reverb Reducer Program	✓
Noise Management	✓
Wind Noise Reduction	✓
Bandwidth Extension	✓
Bandwidth Compression	✓
Intelligent Feedback Preventer	✓
Music Enhancer (3 programs)	✓
Automatic Acclimatization	✓
Smart Automatic Equalizer	✓
Automatic Classifier	✓
XPhone	✓
Remote Services (Basic/portal)	✓
IP68 Rated	✓
ATEX Certified	✓

✓ available - not available



## ATEX CERTIFICATION

Reach hearing aids have an ATEX certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments. They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills.

# REACH B-LI RUGGED

REACH B-LI RUGGED IS OUR TOUGHEST BTE AND CAN HANDLE ALL THE KNOCKS AND SPILLS THAT OCCUR IN EVERYDAY LIFE

Our toughest BTE with Reach technology so the wearer never misses out of the conversation.

RELY ON  
LIFEPROOF  
REXTON



- Rexton App
- Rexton Assist
- TeleCare
- Made for iPhone + Android ASHA
- Headset Mode Tap Control

Rexton's toughest hearing aid, built to withstand harsh conditions, such as rain and showers, and resistant to sweat, salt, and oils:

- > Multi-Voice Focus
- > MyVoice 2.0
- > Bluetooth connectivity via LE Audio, MFi, and ASHA
- > Hands-free calls with compatible iOS and LE Audio smartphones
- > Up to 36 hours of battery life, with up to 32 hours when streaming
- > On-board rocker switch
- > Built-in telecoil lets the wearer hear clearly in noisy environments with loop systems
- > Connexx Charging Station for ease of use with earmolds
- > The earhook can be exchanged for a thin tube, compatible with EarWear 3.0.

#### Unique tough building offering:

- > Water resistant construction
- > Sweat proof
- > Drop proof
- > Scratch resistant
- > IP68 certified
- > ATEX certified

#### Chargeable via



Connexx Charging Station B-M  
Included


#### SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast-charge functionality
- Fixed USB charging cable
- Included power adapter
- LED charging status lights

Colors

Dark Champagne   Granite	Sandy Brown   Granite
Rose Gold   Rose Gold	Black   Black

2 ccm Coupler (All fitting ranges are based on 2 ccm coupler)

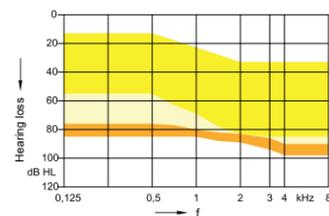

ThinTube 3.0 + EarTip 3.0 Open  
ThinTube 3.0 P + Sleeve 3.0 Powe  
Earhook + Mol

FOG (peak) 60 dB  
OSPL 90 (peak) 133 dB SPL

## REACH B-LI RUGGED

<b>Performance Level</b>	80
Channels/Controls/Programs	48/20/6
<b>Battery</b>	
Lithium-ion battery	✓
Charging time	4 h
Battery runtime (without streaming)	up to 36 h
Battery runtime (incl. 5 h streaming)	up to 32 h
<b>FOG in dB</b>	60
<b>Accessories</b>	
Connexx Charging Station B-M	✓
TV Sound	✓
Smart Mic	✓
Smart Key	✓
<b>App</b>	
Rexton App	✓
<b>EarWear 3.0</b>	
Earhooks (standard & small)	✓
ThinTubes (standard & power)	✓
Domes & Sleeves	✓
Earmold 3.0	✓

## Features

Speech Preservation Technology	✓
Multi-track Processing	✓
Rexton Assist	✓
TeleCare	✓
Motion Sensor	✓
IP68 Rated	✓
Soundpro 2.0	✓
MyVoice 2.0	✓
Direct Streaming	✓
Made for iPhone - Android ASHA- LE Audio	✓
HandsFree/Tap Control*	✓
Auto Volume	✓
Wireless Sync	✓
Auto Stereo iLock	✓
Auto iFocus 360	✓
Noise Reduction	✓
Wind Noise Reduction	✓
Reverb Reducer	✓
Music Enhancer (Live/Recorded/Playing)	✓
T-Coil	✓
XPhone	✓
ATEX Certified	✓

\* iPhone and LE Audio only - not ASHA

## TOUGH DESIGN IN A SUPER-COMPACT SIZE

Reach B-LI Rugged includes an on-board rocker switch, a t-coil, and is compatible with EarWear 3.0. The earhook can be exchanged for a thin tube.

To change the earhook to a thin tube use the tool pictured here. Remove the axle pin and unscrew the earhook until it is completely removed. Place the thin tube and screw it on all the way. The axle pin can be discarded.



## ATEX CERTIFICATION

Reach hearing aids have an ATEX certification which verifies that they are tested and approved to be intrinsically safe for use in potentially explosive environments. They can be safely used in areas where there may be flammable substances in the air such as paint and chemical manufacturing, petrol stations and mills.



# CHARGERS

Rexton offers three different charger types; desktop, Dry & Clean, and portable, so you can choose the one that best suits user needs.



Connexx Charging Station B-M Rugged BTE



Charging Station R RIC

## DESKTOP CHARGING STATION

A highly powerful hearing aid needs a powerful and long lasting charger. Our standard chargers have a galvanic charging system that ensures a lossless charging process.

### SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast-charge functionality
- Fixed USB charging cable
- Included power adapter
- LED charging status lights

### COMPATIBILITY

- Reach Plus BC
- Reach R-Li T
- Reach B-Li Rugged

## DRY & CLEAN STATION (CHARGING+ STATION)

The Charging+Station is our premium charger. Go beyond and offer extra care to the hearing aids with an integrated drying function and UV cleaning light, as it charges.

### SPECIAL FEATURES

- Robust and scratch resistant design
- 30-minute fast-charge functionality
- USB-C power adaptor
- Included power adapter
- Integrated drying functionality
- Integrated UV cleaning functionality

### COMPATIBILITY

- Charging + Station R:
- Reach Plus BC
  - Reach R-Li T

## CHARGERS (CONTINUED)



## TRAVEL CHARGER RIC

This portable charger provides up to 3 charges on the go – without connecting to a power outlet.

### SPECIAL FEATURES

- Holds 3 full charges
- Charges in 4 hours
- 30-minute fast-charge
- Automatic on/off
- USB-C, power adapter
- Integrated Li-ion cell

### COMPATIBILITY

- Reach Plus BC
- Reach R-Li T

## IF CIC CHARGER

Rexton Reach IF CIC Travel Charger gives wearers up to 28 hours of runtime and full charging flexibility with plug-in or wireless charging options.

### SPECIAL FEATURES

- Holds 4 full charges
- USB-C cable
- Wall adapter (included)
- Integrated Li-ion cell
- Qi Charging pad compatible (not included)
- 30-minute fast charging

### COMPATIBILITY

- Reach inoX-CIC Li

# THE COMPREHENSIVE REXTON APP

ONE FOR ALL



The Rexton App allows wearers to control their hearing system conveniently as it turns a smartphone into a discreet remote control for the hearing aids.

Telephone calls are forwarded directly from the smartphone and take priority over other audio signals; streamed music can be listened to in stereo, directly, and without any additional device\*.

Download the Rexton App now!



iOS

Android

## SPECIAL FEATURES\*

- Hearing programs and volumes can be adjusted discreetly
- The volume can be set individually for the right and left hearing aids
- Directionality control
- Remote support
- Real-time battery status
- Make real-time programming changes with Rexton Assist

## DEVICE COMPATIBILITY


Wearers can check if their smartphone is compatible with Rexton hearing aids and app.



Scan the QR code to check smartphone compatibility for:

- Rexton App use
- Streaming phone calls and audio
- Hands-free calling

\*Check the compatibility of hearing aid and smartphone device

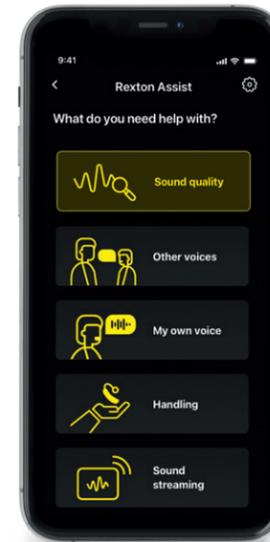


## REXTON ASSIST

Our cloud-based Artificial Intelligence neural network allows wearers to make real-time programming changes when they are on the go. Wearers can personalize their hearing aids' performance when and where challenges are occurring without having to wait for an in-person or virtual appointment.

After receiving solutions curated for that specific situation, wearers provide feedback to the living network so it can evolve based on the wearers' satisfaction level.

Compatible with Reach Plus BC, Reach R-Li T, Reach B-Li Rugged.

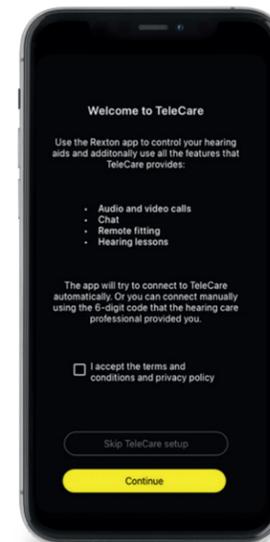

## REXTON TELECARE

With Telecare, Hearing Care Professionals can increase their service journey and provide remote support through:

- Instant messages
- Audio & video calls
- Fine tuning

Compatible with Reach Plus BC, Reach R-Li T, Reach B-Li Rugged.

Reduced Telecare options available to Reach inoX-CIC Li.

# AUDIO STREAMING

Rexton offers audio streaming devices as interfaces to keep users connected for a pleasant hearing experience, no matter what. Compatible with Bluetooth-enabled hearing aids only.



## TV SOUND



The new TV Sound will help the hearing aid user to hear the TV directly in their hearing aids. With an updated look and feel together with a latency less than 30ms for all inputs, the TV sound can be paired with several hearing aid pairs at the same time for a shared enjoyment.

### MODEL FEATURES

- Auto connect to hearing aids already paired
- LED for status showing power, pairing and streaming



## SMART MIC



The Smart Mic enables hands-free telephone calls without the user needing to pick up their smartphone. In challenging hearing environments it can easily be used as an additional microphone, and allows basic hearing aid control functions.

\*Reach R Plus is not compatible with Smart Mic.

### MODEL FEATURES

- Direct streaming from Android smartphones
- Listen to music in stereo quality
- Can be used as a conference microphone

# REMOTE CONTROL

Control the volume of your hearing aid easily and quickly with a touch.



## SMART KEY




With Smart Key, programs and volume can be easily adjusted for the hearing aid. A key lock prevents buttons from being accidentally pressed and a status LED indicates whether the remote control is on or off. Compatible with all Rexton hearing aid models.

### MODEL FEATURES

- Small build to fit on a keychain
- Control via high-frequency sounds

# EARWEAR 3.0

EarWear 3.0 provides a comprehensive solution, which is highly reliable in its acoustic performance and highly durable. Therefore, it improves the cosmetic appeal as well as the physical fit and retention in the ear. So, wearers enjoy the most comfortable, cosmetically appealing solution with consistently great sound.

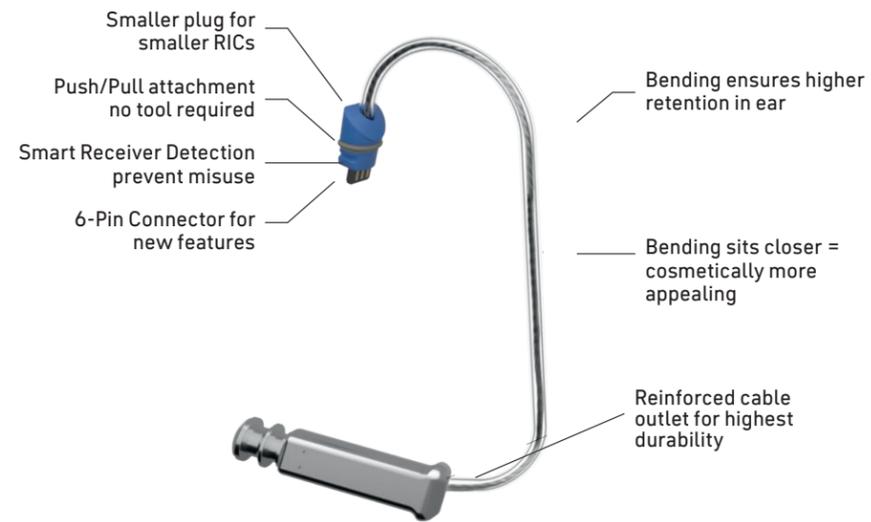

## COMPATIBLE HEARING AIDS

- Reach Plus BC
- Reach R-Li T

## POWER

- S, M, P, HP

## LENGTHS

- 0, 1, 2, 3, 4 and 5

## CONCHA LOCKS

- S, M, P

## EARWEAR EAR TIPS 3.0

Soft tips for a comfortable fit.

### DOMES

Open 4mm, 7mm, 10mm

Closed 8mm, 10mm, 12mm



Tulip 8mm, 12mm

### SLEEVES

Vented XS, S, M, L

Closed XS, S, M, L

Power XS, S, M, L

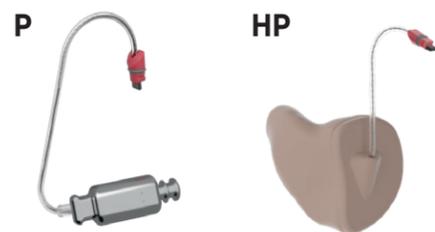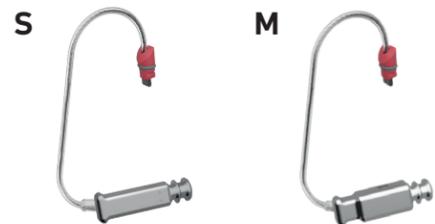

# THINTUBE 3.0

Aesthetically appealing options that are also for severe hearing loss.



## SIZES

- Thin Tube S: 0.9mm for mild to moderate losses
- Thin Tube P: 1.4mm for severe losses

## LENGTHS

- 0,1,2,3 and 4

## COMPATIBLE HEARING AIDS

- B-Li Rugged

miniReceiver 3.0 and ThinTube 3.0 use the same Ear Tips and Sleeves for maximum simplicity and same cable/tube geometry for best cosmetics.




## NANOCARE 3.0 WAXGUARDS

Keep wax from damaging receivers with replaceable waxguards.



## CONCHA LOCK 3.0

Sleeve-style locks are easy to attach and remove.

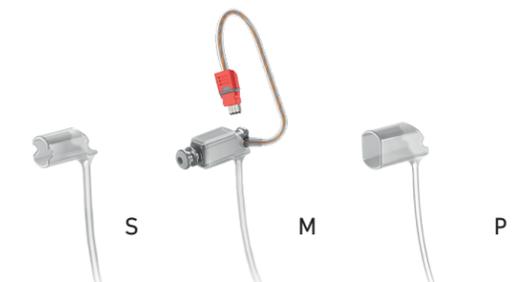



# RELY ON REXTON

It is a tough world out there, and even tougher when your hearing starts to become challenging. At Rexton, we will not let your hearing loss hold you back. Good hearing is critical. From getting the job done, to getting home safely, to being there for the people who count on you, we know what is at stake. That is why we work hard to deliver proven hearing technology in the most reliable way. We understand what reliability means in real life, so we make practical, and durable, easy-to-use products you can count on. Our hearing solutions are Lifeproof, so you can count upon to perform reliably in even the harshest environments. Rexton delivers proven technology in the most reliable way. RELY ON REXTON.

## **Distributor**

WS Audiology, USA Inc.  
101 Wood Ave South, 7th Floor  
Iselin, NJ 08830  
1 (800) 766-4500

[rexton.com](http://rexton.com)

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone or iPad may affect wireless performance. iPad, iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. The Android name, the logo, and other trademarks are property of Google LLC.

The Bluetooth® word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by WSAUD A/S is under license. Other trademarks and trade names are those of their respective owners.

Copyright © 2026 Rexton is a registered trademark of WSAUD A/S.  
10247183 v10 01/26 RX-23704